GESKE
GERMAN BEAUTY TECH

Quickstart Guide

www.geske.com

GESKE Beauty Tech GmbH,
Leipziger Platz 18, 10117 Berlin, Germany
info@geske.com, www.geske.com

SmartAppGuided™
Anti-Aging Skin Renewer | 9 in 1

Short Introduction


Scan for the full instructions in your language

GESKE German Beauty Tech App
Available in 45 Languages

Apple and the Apple logo are trademarks of Apple Inc.
Google Play and the Google Play logo are trademarks of Google LLC.

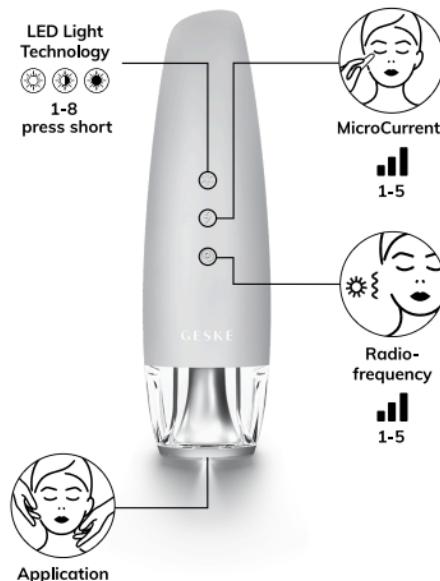EN

Read the safety information before first use. Keep both instructions for use.

YOUR PRODUCT AT A GLANCE

Your product is a multifunctional beauty device with radio frequency, microcurrent, and LED technology. The two most intensive technologies, radio frequency and microcurrent, can be safely used up to 3 times per week. The LED technology can be used daily. For maximum safety and optimal results, always use your personalized skincare routine, which was created by the GESKE app.

Download the GESKE app and get your personalized routine.

SPECIFICATIONS:

Material: ABS, acrylic, PC, silicone, stainless steel | **Storage Temperature:** -5 – +60 °C / 23 – 140 °F | **Battery:** Li-Ion 3.7V | **Battery Capacity:** 1,000 mAh | **Maximum Charging Current:** 730 mA | **Charging Voltage:** 5V | **Charging Time:** approx. 3 hours | **Operating Time:** Operating time differs depending on activated technologies & intensities | **Interface:** USB | **Programs:** LED Decorative Light: 1 color (White), Vibration: 5 speed levels, RF: 5 intensities, MicroCurrent: 5 intensities, LED Light Technology: 8 colors

IMPORTANT SAFETY INFORMATION! READ BEFORE USE!

- Information on radiofrequency:
- This product uses multipolar radiofrequency, which targets deep layers of the skin and superficial layers of fatty tissues in a non-ablative way, e.g. without damaging the epidermis
- To reduce the risk of potential side-effects and to avoid over-heating or burning of the skin, always keep the device in motion and do not use more frequently than recommended above
- Always start with the lowest intensity and closely monitor your skin's reactions before increasing the intensity level
- Several sessions at a lower intensity provide more substantial effects than a single session with higher intensity
- Consistent use of the device is required to see continuous improvements and results
- Expected side effects:
 - During or shortly after usage you may experience erythema (redness of the skin) and edema (swelling of the skin). Those side-effects are harmless and should disappear within a few hours, at maximum after 24–48 hours after usage.
 - After using the device on your face you may experience headaches which will disappear within 24–48 hours after usage.
 - Unexpected side effects:

• Overtreatment with RF (too often or too intense) may cause depressions on the skin, burnings, scarring, and damages on the fatty tissue or epidermis. To avoid this, only use the device as recommended

• After intensively using the device on your neck, the skin below your ear may feel numb. This rare side effect is temporary, but may take some weeks to completely subside

• Further safety information:

• Do not use the product on irritated or injured skin or open wounds (after any kind of injury), in the case of bacterial or viral infections of the skin (either in the usage area and elsewhere)

• This product is not intended to treat medical skin conditions or diseases. If you have medical skin conditions or diseases, consult your dermatologist before using this product

• Do not use the product if it is 1) untypically warm, 2) mechanically damaged or deformed, 3) discolored. Due to RF technology, it is expected that the device becomes warm during usage

• If pain and/or discomfort exceeding the expected side-effects occurs during use, discontinue use and seek medical advice

• The eye area is very sensitive: Be extra careful in this area and avoid any contact of the product with the eyelids or the eye itself

• Only use the product as recommended in your personalized routine

• For hygiene reasons, we recommend you do not share your product with anyone else

• The product is intended for cosmetic use on the face. GESKE assumes no liability for harmful consequences arising from the use on other body areas or other improper use, connection to unsuitable voltage sources, contaminated conducting liquid or balls

• Always make sure the plug and socket are completely dry before charging. Otherwise, electric shocks, a short circuit or fire may result. The product must not be used while connected to the power outlet

• The product can emit intense light, avoid applying it in the direction of or near an open eye

The product is not suitable for use under the following conditions:

- If you have an implanted medical device, metal or other electronic device such as a pacemaker
- After laser treatment or chemical peeling
- After facial surgery, both cosmetically and medically indicated
- Do not use if you are known to be hypersensitive to heat (e.g. heat urticaria)
- If a feeling of (too) intense heat occurs during usage, it must be stopped immediately

In the following cases, a consultation with a physician is necessary as to whether the product is suitable for you to use:

- If you have a heart disease
- If you suffer from an autoimmune disease or epilepsy, consult a specialist
- In case of pregnancy. The safety of use during pregnancy has not been established
- In case of anamnestically known, clinically healed Herpes Simplex on the face: In rare cases, the use of an intensive cosmetic device can reactivate a pre-existing Herpes Simplex Infection of the face (especially in the area of the lips)

Further information on use:

- Before using the MicroCurrent Pro Face-Lift Technology and the RF Heated Wavelengths, the product must be completely dry
- Use of the MicroCurrent Pro Face-Lift Technology and the Advanced RF Heated Wavelengths only in combination with the GESKE Hydrating MicroCurrent Gel
- Do not use the MicroCurrent Pro Face-Lift Technology and the Advanced RF Heated Wavelengths on the chest area, around the eyes, throat, genital/groin area, or on red and inflamed skin, over birthmarks and enlarged blood vessels
- Nothing is known about the long-term effects of MicroCurrent and RF
- As a result of stimulation of the optic nerve, a flashing light may be perceived during use. If this persists when not using the product, seek medical advice
- A slight tingling sensation may occur during application. Decrease the intensity of the MicroCurrent or the RF to reduce or eliminate tingling
- Electronic monitoring devices may not function properly while the device is in use
- Never use the product while driving or operating heavy machinery
- Do not poke any sharp objects into the product's openings
- The product may only be used with a safety extra-low voltage (SELV) adapter
- We recommend using standard IEC60335 power supplies to charge the device
- The safety and effectiveness of the product has never been assessed in conjunction with non-GESKE skincare products. We recommend using the product only with our microcurrent products. Use with other products at your own risk
- The recommended frequency of use of the device depends on the technology selected. Radiofrequency should not be used more than three times a week. The other technologies can be used safely up to a maximum of once a day
- You can also use the different technologies one after the other. The Advanced RF Heated Wavelengths should always be used last to avoid skin irritation or injury. Please note that applying several technologies in succession increases the risk of skin irritation such as redness or slight swelling
- For 24 hours after an intensive session with MicroCurrent Pro Face-Lift Technology or Advanced RF Heated Wavelengths, you should avoid cosmetic treatments with other GESKE devices as well as sports and sauna visits. Use sunscreen with a high SPF in the

days following the session

- Keep product out of reach of children

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